



*Technology is the art of creating usable
realities from simple thoughts.*

T. F. Lehnert



Welcome to **Baumann-Dental**

We are delighted to present to you our new catalogue full of great new ideas and useful products! We've placed the priority on clarity, understanding and presentation, in order to give you clear expectations when it comes to quality.

Further information, brochures and instructions on many products can also be found in the service area of our website:
www.baumann-dental.de

Of course, nothing can beat personal advice. Our expert team will be happy to take the time to respond to your questions and provide you with detailed information about all our products.

Fon +49 7236-9 33 69 0

We look forward to your call!

Innovation and precision

Baumann-Dental has been operating successfully on the market for almost 30 years. What started very small in 1988 has grown over the years to become an internationally recognised company in the sector.

We were able time and again to come up with developments that have since been distributed throughout the world. Our Adesso Split magnetic plate system, for example, is still the gold standard when it comes to precision when coordinating articulators, even after over 25 years.

Taking a good look and simply making products a bit better than others is what Baumann Dental does: precise and innovative. This year, too, you can find a range of new

developments in our catalogue. We would like to thank you, our customers throughout the world, who always give us constructive feedback and help us to further develop our products and continuously optimise them.

Always be sure to use the original from Baumann Dental!

Best regards,
The Management Board at Baumann Dental

Karsten Baumann
Gerhard Buchenroth
Wolfgang Dawedeit



Gerhard Buchenroth

Karsten Baumann

Wolfgang Dawedeit

Important developments at Baumann-Dental GmbH

1988

Company founding

Baumann Dental GmbH was founded by Karsten Baumann in Mühlacker/Enzkreis at the start of 1988.

1989

adesso Split® magnetic plate system

The gold standard, both then and now, when it comes to precision and the coordination of articulators.

1998

Model System 2000®

This model system is characterised by its high precision, speed and ease of use.

2005

Modelsystem comBi flex®

The affordable alternative to the Zeiser⁴/ Giroform¹ system. Take a good look!

2009

adesso MULTISPLIT® magnetic plate system

The cleverest idea since the development of magnetic plate systems

2011

ArTO® articulators

Want robust and uncomplicated instruments? Then Arto is perfect for you!

2015

CoNECT mounting plates

Get to know a new kind of model mounting: absolutely reliable and very economical.

01 We work for you

Welcome to Baumann-Dental GmbH	Page 5
Important developments at Baumann Dental GmbH	Page 7
Table of contents	Page 8-9

02 Articulators

Theme image articulators	Page 10-11
Arto articulators	Page 12-17
Articulator accessories & replacement parts	Page 18-19
Articulator extensions	Page 20-21
Mounting articulators	Page 22-23
Physio-Logic articulator	Page 24-27

03 Coordination

Theme image Adesso Split	Page 28-29
Adesso Split magnetic plates	Page 30-32
Adesso Split milling base	Page 33
Adesso Split accessories	Page 34-37
Adesso Split calibration key	Page 38-41
Theme image Adesso Multisplit	Page 42-43
Adesso Multisplit magnetic plates	Page 44-46
Adesso Multisplit milling base	Page 47
Splitex ¹ compatible magnetic plate system	Page 48-49
KaVo ³ Split compatible magnetic plate & calibration key	Page 50-51

04 Model mounting

Theme image model mounting	Page 52-53
Connect mounting plates for Adesso Split & Splitex ¹	Page 54-55
Mounting plates for Adesso Split, Splitex ¹ & KaVo ³	Page 56-57
Premium mounting plates for Adesso Split & adhesive disks	Page 59

05 Model making

Theme image model manufacture	Page 60-61
Model System 2000	Page 62-63
Model System 2000 starter kit	Page 64-65
Model System 2000 components & accessories	Page 66-69
Duet plate system	Page 70-71
Duet plate system starter set & plates	Page 72
Model System 2000 Evo starter kit	Page 74-75
Model System 2000 Evo components & accessories	Page 76-77
Combiflex model system Premium	Page 78-79
Combiflex model system Basic & Splitcast	Page 80
Accessories for Combiflex, Zeiser ⁴ & Giroform ¹	Page 81-83
Zeiser base plates ⁴ & accessories	Page 85
Quadra Base model system	Page 86-87
Model base moulds & matrixes	Page 88
Xpress pin system & accessories	Page 89
Occlusion spray / silicone spray / membrane boxes	Page 90
Splitcast items	Page 91

06 Grinding Equipment

Grinding belts for dry trimmers	Page 92
Mandrels & grinding belts / tungsten carbide cutter	Page 93

07 Duplicating Technology

Duplicator	Page 94
Basic duplicating cuvette	Page 95

08 Transfer equipment

Transfer stands / transfer carriages	Page 96-97
Transfer table / bite fork support / accessories	Page 98

09 Workshop service & special conversions

Workshop service & special conversions	Page 99
--	---------

10 Sales & service / T&Cs

Sales & service / Legal information	Page 100
General Terms and Conditions of Baumann-Dental GmbH	Page 101



AN ARTO FOR ALL EVENTUALITIES

The Arto articulator range offers impressive performance at a low price. With a consistent classic fully aluminium construction, these non-arcon articulators fulfil all your requirements in terms of stability and precision in every-day lab work.

Arto articulators are available in the construction height 116 mm (construction height of the 'old' Artex¹), as well as in the construction height 126 mm (Artex¹ carbon articulator). Depending on the requirement, three different versions of the Arto concept are available to you. From the classic average value articulator Arto 1 to the partially adjustable Arto 3, you are sure to find the right equipment.

The articulators are supplied from the factory with model plates made from plastic. To make the most of the full flexibility of the articulators, we recommend incorporating an Adesso Split or Adesso Multisplit magnetic plate system. Precise, stable, economical... Arto from Baumann Dental!

ArTO 1® The beginner

A classic average value articulator, robust and long-lasting, with the high precision you expect.

The features in detail:

- Available in the construction heights 116 or 126 mm
- Bennett angle (fixed) 12°
- Condylar inclination (fixed) 35°
- Table support 45° (removable)
- Plastic mounting plate (1 pair premounted)

Suitable for the mounting of magnetic plate systems:

- Adesso Split Set / Artek¹ / Arto / #10000
- Adesso Multisplit Set / Artek¹ / Arto / #10200
- Splitex¹ compatible magnetic plate set / Artek¹ / Arto / #10010

12400 Arto 1, 116 mm / incl. 2 pcs. plastic model plates

12401 Arto 1, 126 mm / incl. 2 pcs. plastic model plates



Tip: As all Arto articulators are already fitted with plastic mounting plates, you can begin work immediately after unpacking them. However, the articulators only unfurl their full potential when assembled with our magnetic plate systems, i.e. the tried-and-tested Adesso Split or the exceptionally flexible Adesso Multisplit.



Figure shows: *Arto 1* tilted with a construction height of 116 mm



Figure shows: *Arto 1* with a construction height of 116 mm



Figure shows: **Arto 2** with a construction height of 116 mm

ArTO 2® The climber

Our best-seller! Fully sufficient in terms of its functions for most dentistry work.

The features in detail:

- Available in the construction heights 116 or 126 mm
- Bennett angle (adjustable) 0° to 20°
- Condylar inclination (fixed) 35°
- Table support 45° (removable)
- Plastic mounting plate (1 pair premounted)

Suitable for the mounting of magnetic plate systems:

- Adesso Split Set - Artex¹ / Arto / #10000
- Adesso Multisplit Set - Artex¹ / Arto / #10200
- Splitex¹ compatible magnetic plate set / Artex¹ / Arto / #10010

12500	Arto 2 - 116 mm / incl. 2 pcs. plastic model plates
12501	Arto 2 - 126 mm / incl. 2 pcs. plastic model plates



Arto 2 Bennett angle adjustable from 0°-20°. Non-slip knurled screws enable secure fixation.



Figure shows: **Arto 2** with a construction height of 126 mm

! *Tip: Do you still have usable Artex¹ articulators of an older type? The anterior guidance pins, holders and plates are often very worn in these articulators and can no longer be used. Save a lot of money and replace these defective parts with our Arto replacement parts, as some of these are compatible. Ask us!*

ArTO 3® The pro

The top model in the Arto range, just as precise and robust as its little brothers, but comparatively good value for money.

The features in detail:

- Available in the construction heights 116 or 126 mm
- Bennett angle adjustable 0° to 20°
- Condylar inclination adjustable -15° to +60°
- Arbitrary pins for direct face bow adaptation
- Table support 45° (removable)
- Plastic mounting plate (1 pair premounted)

Suitable for the mounting of magnetic plate systems:

- Adesso Split Set - Artex¹ / Arto / #10000
- Adesso Multisplit Set - Artex¹ / Arto / #10200
- Splitex¹ compatible magnetic plate set - Artex¹ / Arto / #10010

12600 Arto 3, 116 mm / incl. 2 pcs. plastic model plates

12601 Arto 3, 126 mm / incl. 2 pcs. plastic model plates

! *Tip: With direct adaptation of a face bow to the arbitrary pins of the articulator, the condylar inclination should be set to 60 degrees so that the registration is in the correct position. For more comfortable working, we recommend our transfer accessories on pages 96-98*

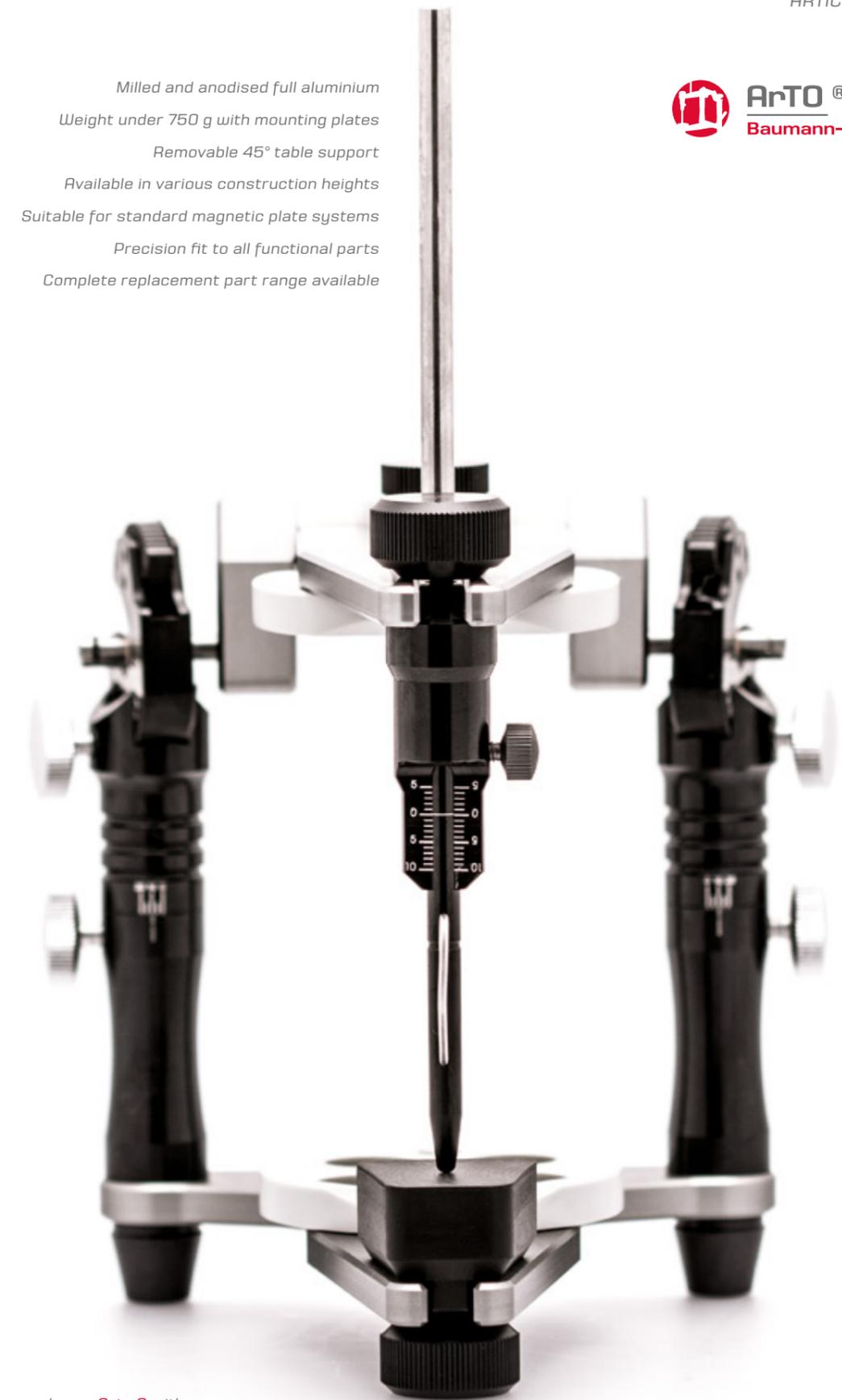


Figure shows: Arto 3 with a construction height of 116 mm



Arto 3 Condylar inclination adjustable from -15° to +60°. Non-slip knurled screws enable secure fixation




Arto 3 Bennett angle adjustable from 0° to 20°



Arto 3 Arbitrary pins for direct face bow adaptation



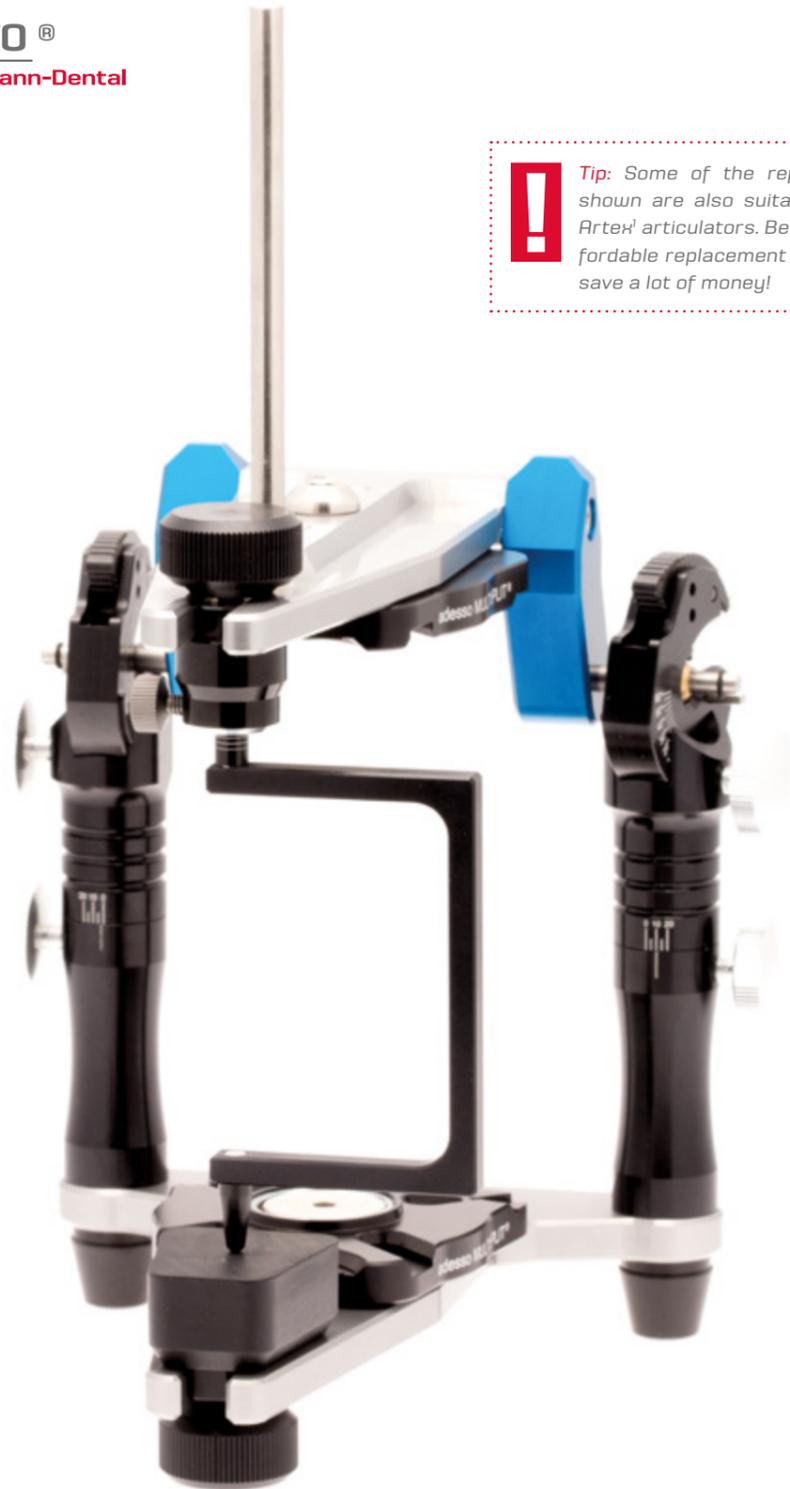

! *Tip: Some of the replacement parts shown are also suitable for use with Artex¹ articulators. Benefit from our affordable replacement parts prices and save a lot of money!*

Figure shows **Arto 3** in the construction height 126 mm with **free view anterior guidance pin**, for a better overview in prosthetics

ArTO® & Artex¹ Replacement parts & accessories

Here you can find the most important replacement parts for our Arto articulators. Naturally, you can also get all other replacement parts as well - just ask us!

12003	Incisal needle for articulators / without illustration
12110	45° table support for Arto / without illustration
12111	Stabilising pin for Arto / without illustration
12112	Anterior guidance pin for Arto / 116 mm
12112-02	Anterior guidance pin incl. holder for Arto / 116 mm
12113	Anterior guidance pin for Arto / 126 mm
12113-02	Anterior guidance pin incl. holder for Arto / 126 mm
12151	Free view anterior guidance pin for Arto 116 mm / Artex ¹
12152	Free view anterior guidance pin for Arto 126 mm / Artex ¹ Carbon
12305	Anterior guidance table for Arto 0° / incl. screw
12150	Anterior guidance pin for Artex ¹ / 116 mm
15800-01	Anterior guidance pin for Artex ¹ Carbon / 126 mm



ArTO® model plate set

Can be used with Arto, Artist¹ and Artex¹ articulators (not suitable for Artex¹ Carbon)

12000	Model plate set for Arto / Artex ¹ / 2 pcs.
12001	Model plate screw set for Arto / 2 pcs.
12002	Model plate screw set for Artex ¹ / 2 pcs.



*...do you still have questions?
We look forward to your call...*
Fon +49 7236-9 33 69 0

Articulator extensions

Articulator extension Artex¹ to Artex¹ Carbon for the models AP / AM / AR / AL / AN (ARCON) & TS / TR (NON ARCON)

This conversion kit makes it possible to raise older Artex¹ arcon models to the new construction height of the Artex¹ Carbon articulator without influencing symmetry. Suitable for Artex¹ AP, AL, AR, AN and similar models. Also for use with Artex¹ non-arcon models with immediate side shift function (TS / TR).

15810 Articulator heightening modification for Artex¹ Arcon / AM, AP, AL, AR, AN, TS, TR / 116 to 126 mm



Figure shows articular extension #15810 used with Artex¹ TS / TR



Figure shows articular extension #15810 used with Artex¹ arcon model

Articulator extension Artex¹ to Artex¹ Carbon for the models NK / N / TK and similar

This clever heightening modification enables you to bring your old Artex¹ to the construction height of the new Carbon articulator in just a few minutes. Suitable for self-assembly.

15800 Articulator extension for Artex¹ / NK, N, TK / 116 to 126 mm



Figure shows: Parts supplied for extension #15800 for Artex¹

Other articulator extensions

Almost all available heightening modifications are suitable for self-assembly. For Panadent articulators, we offer a conversion service with various options - please just ask us!

15000	Extension for SAM ² 2 to SAM ² P / 105 to 121 mm
15100	Extension for Balance / 95 to 116 mm
15200	Extension for Balance / 105 to 116 mm
15300	Extension for Dentatus / 104 to 116 mm
15400	Extension for Dentatus / 104 to 121 mm
15500	Extension for Panadent / various options, on request!
15600	Extension for Rational / 105 to 116 mm
15700	Extension for SAM ² 1 to SAM ² P 'old' / 106 to 121 mm
15701	Extension for SAM ² I to SAM ² P 'new' 106 to 121 mm
15900	Extension for Artex ¹ Carbon CN, CT, CP / 126 to 136 mm. Suitable for integration in Artex ¹ Carbon non-arcon models to create more space for model mounting with the Hip-Mount [®] device.

...do you still have questions?
We look forward to your call...
Fon +49 7236-9 33 69 0

Mounting articulator

incl. magnetic plate system

adesso MULTISPLIT® or adesso split®

Our plaster equipment has been constructed uncompromisingly for tough daily lab use and can still be used reliably for model mounting, even after many years. Just a few moving parts, a very precisely made axis and the special coating help to make the equipment extremely long-lasting. Delivery includes Adesso Split or Adesso Multisplit magnetic plate system from Baumann-Dental.

12100-01	Mounting articulator for Arto 116 mm / Artek¹ incl. Adesso Split
12100-02	Mounting articulator for Arto 126 mm / Artek¹ Carbon - incl. Adesso Split
12100-07	Mounting articulator for SAM² 2P/3 - incl. Adesso Split
12200-01	Mounting articulator for Arto 116 mm / Artek¹ - incl. Adesso Multisplit
12200-02	Mounting articulator for Arto 126 mm / Artek¹ Carbon incl. Adesso Multisplit
12200-07	Mounting articulator for SAM 2P/3 incl. Adesso Multisplit
12010	Plaster protection plate / 1 pc. /

All Mounting articulators are supplied incl. 2 plaster protection plates #12010.




When using Adesso Split mounting plates on the Adesso Multisplit, we recommend the use of the plaster protection plates supplied to prevent the valuable magnetic plates from being soiled.



The integrated 3-point support enables secure and reliable overhead mounting of the models

...do you still have questions?
We look forward to your call...
Fon +49 7236-9 33 69 0

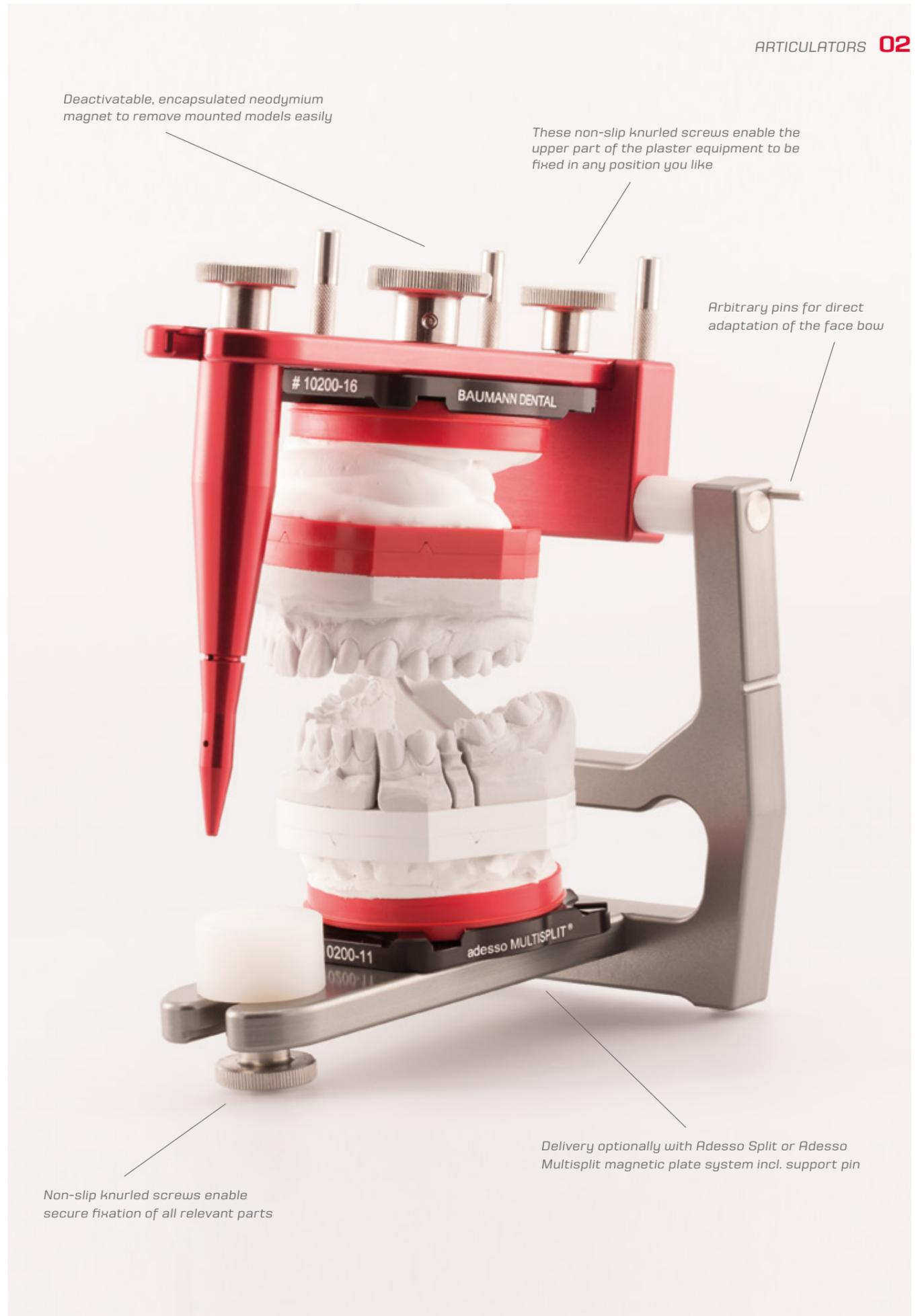

Deactivatable, encapsulated neodymium magnet to remove mounted models easily

These non-slip knurled screws enable the upper part of the plaster equipment to be fixed in any position you like

Arbitrary pins for direct adaptation of the face bow

Non-slip knurled screws enable secure fixation of all relevant parts

Delivery optionally with Adesso Split or Adesso Multisplit magnetic plate system incl. support pin

Physio-Logic articulator

A completely new articulator design



Institut für Temporo-
Mandibuläre Regulation
www.itmr.info

Rainer Schöttl
D.D.S (USA)

*Physio-Logic
Articulation®*
by Rainer Schöttl D.D.S (USA)

Rainer Schöttl

The design gives it away. To construct this myocentric articulator, we had to leave conventional paths far behind...

A completely new concept was developed that is characterised by its logical simplicity and that does away with extensive setting options.

The optimised construction height enables the simple mounting of the models with the Hip-Mount® (www.hip-mount.de), as well as correct positioning via a vector analyser.

For the first time, combined movements, as performed by the human lower jaw, can be specifically replicated in an articulator. In contrast to a conventional articulator, in which a condyle fundamentally stops at its stop point during a laterotrusion movement, the Physio-Logic articulator offers a functionally reversed condyle movement thanks to its special construction in order to set off the greatest possible occlusal compass. The result is maximum trouble-free occlusion.

12603	Physio-Logic articulator / without magnetic plate set
10200	Adesso Multisplit Set - Artex ¹ / Arto / Physio-Logic
10265	Adesso Multisplit spacer plate 10 mm
10266	Adesso Multisplit spacer plate 20 mm

...do you still have questions?
We look forward to your call...
Fon +49 7236-9 33 69 0

Figure shows the Physio-Logic articulator #12603 used with Hip-Mount®. (The Hip-Mount® and magnetic plates shown are special equipment and are not supplied in the scope of delivery.)



Figure shows the Physio-Logic articulator #12603 with Adesso Multisplit System #10200 and 20 mm spacer plate #10266. Prevent an unnecessarily high plaster base by using the spacer plates.

Physio-Logic
Articulation®
by Rainer Schöttel D.D.S (USA)


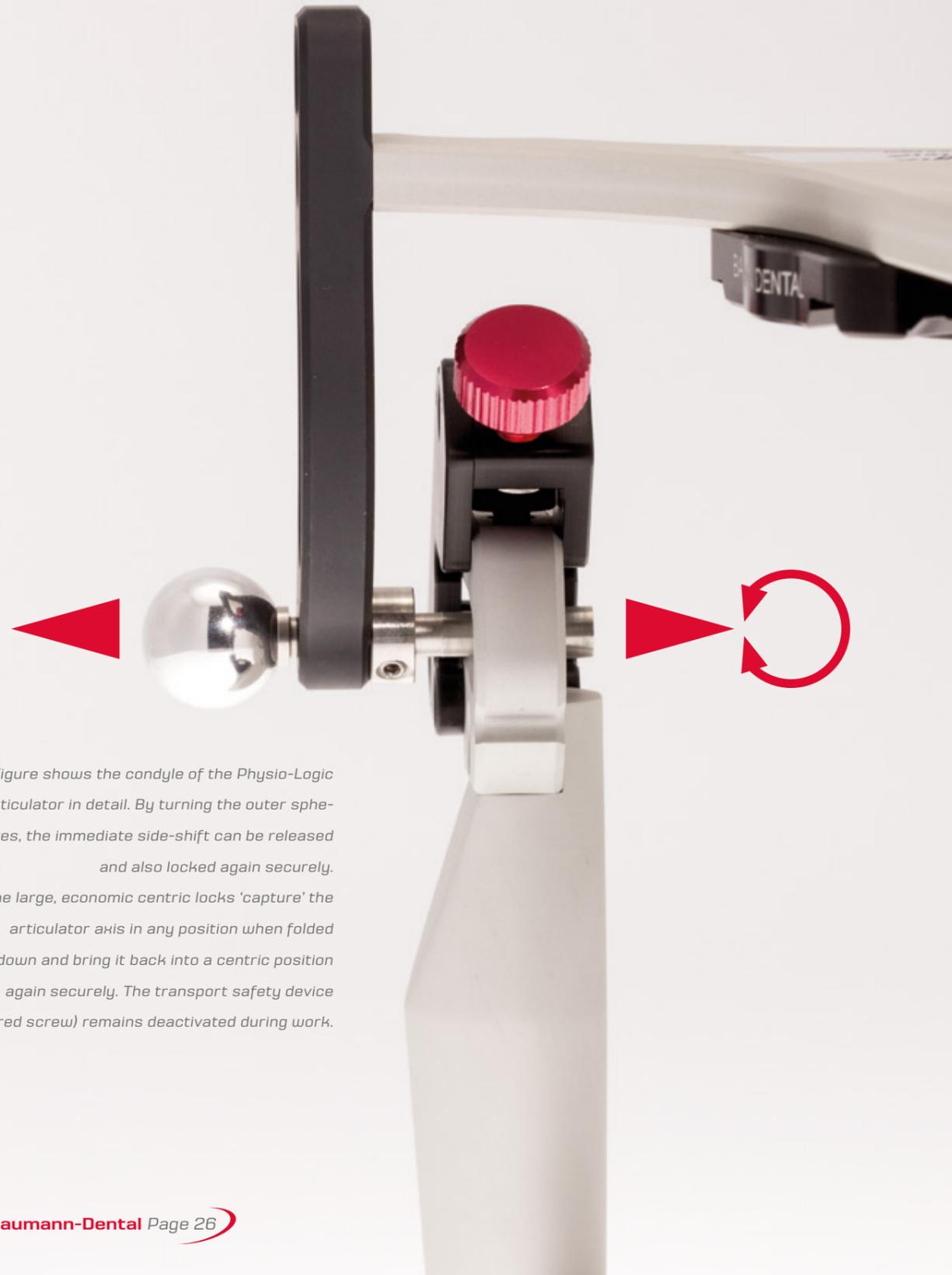



Figure shows the condyle of the Physio-Logic articulator in detail. By turning the outer spheres, the immediate side-shift can be released and also locked again securely. The large, economic centric locks 'capture' the articulator axis in any position when folded down and bring it back into a centric position again securely. The transport safety device (red screw) remains deactivated during work.

Physio-Logic condyle function

The condyles of the Physio-Logic articulator work in a completely different way to the condyles of a classic articulator. The control system is greatly simplified and is designed based on the scanning of a maximum occlusal compass. Any interfering contacts can thus be localised easily. Thanks to the new, reversed path of the condyle spheres, there is no need for an adjustable Bennett angle; this is a fixed part of the movement profile. The easy-to-use immediate side-shift offers an additional option.

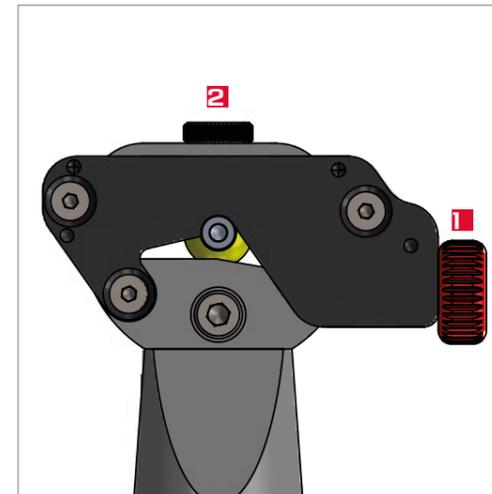


Condyle in a completely locked position with activated transport safety device **1** and centric screw **2**

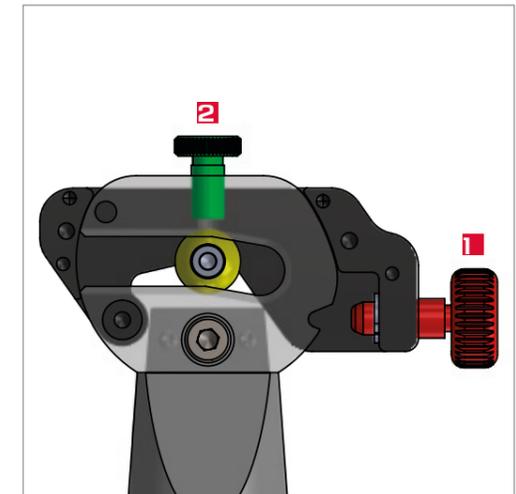

Condyle with deactivated transport safety device **1** and centric screw **2**



Condyle completely opened, condyle sphere is free-moving, centric screw is deactivated.

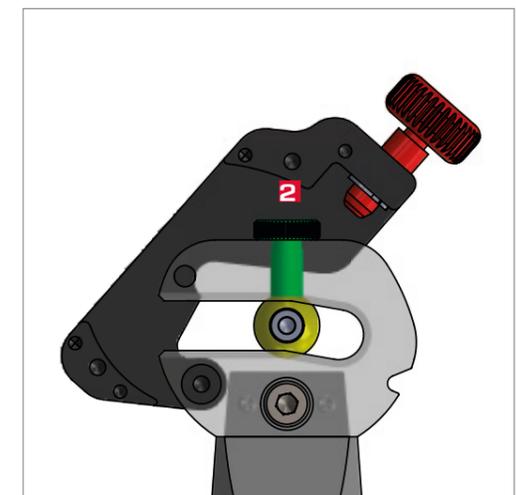

Condyle opened and centric screw **2** activated. Condyle sphere fixed in centric position. Protrusion and retrusion are locked, only movements away from the median level are possible (laterotrusion side).

adesso split® The original

The classic amongst magnetic plate systems

Today, just like 25 years ago, the Adesso Split system sets standard in terms of precision and repeat accuracy. The sophisticated design creates the contact to the mounting plate exclusively via the conically arranged contact surfaces. The mounting plate centres itself independently each time it is fixed in place. The Adesso Split system offers

consistent coordination of various articulator types and thus saves time and money during everyday lab work. The models can be easily removed and used in another, coordinated articulator with no loss of precision. The great variety of available symmetries and important accessories also offers the right solution for your articulators.



Only the conical guiding elements create contact between the Adesso Split base and mounting plate; the model base centres itself perfectly. Soiling on the mounting plate does not impair the fit.

*...do you still have questions?
We look forward to your call...
Fon +49 7236-9 33 69 0*





Figure shows products supplied
Adesso Split Set Artex¹ Carbon #10000-C

adesso split[®] base plate

Mounting takes place in the upper part of the articulator; the plate is screwed in place. Supplied without magnet #10020.

10000-10	Adesso Split base plate / Artex ¹ / Arto
10000-C-10	Adesso Split base plate / Artex ¹ Carbon
10001-10	Adesso Split base plate / Balance
10002-10	Adesso Split base plate / Denar
10003-10	Adesso Split base plate / Dentatus ⁷
10004-10	Adesso Split base plate / KaVo Protar ³
10005-10	Adesso Split base plate / Panadent
10007-10	Adesso Split base plate / SAM ²
10008-10	Adesso Split base plate / Stratos
10009-10	Adesso Split base plate / KaVo Protar ³ Evo
10015-10	Adesso Split base plate / SAM ² = Artex ¹
10017-10	Adesso Split base plate / Reference SL
10064	Adesso Split base plate +10 mm / Artex ¹ Carbon 126 to 116 mm
101010-10	Adesso Split base plate / Dentalhoby / without magnet
10020	Adesso Split magnet / without illustration / graduated prices! for use with Adesso Split, Adesso Multisplit plates and other magnetic plate systems, e.g. Splitex ¹ or Quicksplit ⁶

adesso split[®] systems

All available systems are supplied including magnetic plates #10020. Exception: the systems for the Carbon articulators from AmannGirrbach do not contain any magnets, as these are already supplied with the articulator.

10000	Adesso Split Set / Artex ¹ , Arto
10000-C	Adesso Split Set / Artex ¹ , Carbon / 126 mm / without magnets
10000-CD	Adesso Split Set / Artex ¹ Carbon / 116 mm / without magnets / brings Artex ¹ Carbon to the 'old' Artex ¹ construction height 116 mm
10001	Adesso Split Set / Balance
10002	Adesso Split Set / Denar
10003	Adesso Split Set / Dentatus ⁷
10004	Adesso Split Set / KaVo Protar ³
10005	Adesso Split Set / Panadent
10007	Adesso Split Set / SAM ²
10008	Adesso Split Set / Stratos
10009	Adesso Split Set / KaVo Protar ³ Evo
10015	Adesso Split Set / SAM ² = Artex ¹ / brings SAM ² to the symmetry Artex ¹ construction height 116 mm
10017	Adesso Split Set / Reference SL
101010	Adesso Split Set / Dentalhoby

! *Tip for your order:* our magnetic plate systems are generally designed for the calibration of articulators. You always receive a base plate with thread (for screwing to the upper part of the articulator), as well as an adjusting base plate without thread (for fixing in place with adhesive and calibration key in the lower part of the articulator). If you do not want to calibrate and just want to use the superior splitcast properties of Adesso Split or our other magnetic plate systems, simply use a base plate for screwing at the top and bottom of the articulator. All base plates that can be screwed in place have the identifier '-10' at the end of the item number. As individual base plates are always supplied without magnets, you also need a magnetic plate #10020 for each one. Please note on ordering.

...do you still have questions?
We look forward to your call...
Fon +49 7236-9 33 69 0



Figure shows: *Adesso Split adjusting base plate #10000-11 in the Arto articulator*

adesso split® adjusting base plate

Mounting takes place in the lower part of the articulator; the plate is adhered in place. Supplied without magnet #10020.

10000-11	Adesso Split adjusting base plate / Artex ¹ , Arto
10000-C-11	Adesso Split adjusting base plate / Artex ¹ , Carbon
10001-11	Adesso Split adjusting base plate / Balance
10002-11	Adesso Split adjusting base plate / Denar
10003-11	Adesso Split adjusting base plate / Dentatus ⁷
10005-11	Adesso Split adjusting base plate / Panadent
10007-11	Adesso Split adjusting base plate / SAM ²
10008-11	Adesso Split adjusting base plate / Stratos
10009-11	Adesso Split adjusting base plate / KaVo Protar ³ / Protar ³ Evo
10015-11	Adesso Split adjusting base plate / SAM ² = Artex ¹
10017-11	Adesso Split adjusting base plate / Reference SL Poss. required adapter plate #10017-02
101010-11	Adesso Split adjusting base plate / Dentalhoby
10020	Adesso Split magnet / without illustration / graduated price! for use with Adesso Split, Adesso Multisplit plates and other magnetic plate systems, e.g. SpliteX ¹ or Quicksplit ⁶

adesso split® milling base & milling base plate

The newly designed milling base deserves your attention: Particular value has been placed on an extremely flat construction in order to process higher models including base. You can determine the working height yourself beforehand. The milling base is exceptionally stable and easy to use.

The features in detail:

- Extremely flat construction
- Can be rotated in several directions by up to 90°
- Adjustable height
- Magnetic base plate
- Also for Adesso Multisplit and Model System 2000

10039-01	Adesso Split milling base
10041	Adesso Split milling plate with extra wide edge for more security when clamping / incl. magnet



Figure shows the *milling base* with Adesso Split base plate and the milling plate



adesso split® accessories

Important accessories for magnetic plate systems.

10021	Adesso Split magnet lifter for damage-free removal and application of the magnets #10020
11000	Adesso Split® synchronize adhesive / 2-component adhesive to affix the adjusting base plate in the lower part of the articulator. The special recipe guarantees the greatest possible freedom from shrinkage and excellent resistance to ageing.
11001_M6	Adesso Split special screw for SAM ¹ / KaVo ³
11001-OK	Adesso Split special screw for Artex ¹ / Arto / upper jaw
11001-UK	Adesso Split special screw for lower jaw Artex ¹ & Physio Logic
11006	Adesso Split locking device for SAM ² protects the screw connection to articulators against access by unauthorised users. Parts supplied: Screw locking device and Allen key
11006 / Balance	Adesso Split locking device for Balance protects the screw connection to articulators against access by unauthorised users. Parts supplied: Screw locking device and Allen key
10020	Adesso Split magnet / graduated prices! for use with Adesso Split, Adesso Multisplit plates and other magnetic plate systems, e.g. Splitex ¹ or Quicksplit ⁶

adesso split® accessories

Special plates and accessories for Adesso Split and other accessories can be found in the 'Transfer equipment' section

10043	Adesso Split positive plate / suitable for subsequent conversion of a transfer table
11008	Adesso Split set up template carrier for Artex ¹ / Arto 116 mm
11008-C	adesso split set up template carrier for Artex ¹ Carbon / Arto 126 mm (without illustration)

adesso split® spacer plate 10 mm

Magnetic plate #10020 supplied.

10065	Adesso Split spacer plate 10 mm / Artex ¹ Carbon
10064	Adesso Split base plate +10 mm / Artex ¹ Carbon / without magnet to reduce the Carbon articulator to the construction height 116 mm



Figure shows the Arto3 articulator with Adesso Split system #10000 and Adesso Split 10 mm spacer plate #10065.



adesso split® master plate / matrix

To produce secondary plates from plaster or plastic. To maintain the matrices and for insulation, we recommend our silicone care spray #11009.

10040	Adesso Split master plate (single mould) incl. matrix
10045	Adesso Split master plate (triple mould) incl. matrices
11002	Adesso Split matrix / red
11003	Adesso Split matrix 15° raked
11009	Adesso Split silicone care and insulation spray / 400 ml



Our silicone spray is perfectly suited to caring for your valuable articulators, magnetic plate systems, etc. It is also suitable for insulating your Adesso Split matrices.



! Tip: Never allow matrix surface to come into contact with plaster solvent as the surface is destroyed by plaster solvent!

...do you still have questions?
We look forward to your call...
Fon +49 7236-9 33 69 0



CALIBRATION KEY

Perfect working conditions in your articulator with adesso split®

The Adesso Split calibration key is an adjustment and testing tool in one. Thanks to the patented construction, the key can be divided in two for testing purposes. In contrast to one-piece calibration keys, the Adesso Split calibration keys enable even the smallest inaccuracies to be identified quickly and easily. Even defective condyle housings or distorted frames on your articulators can be identified quickly using the Adesso Split calibration key.

When performing tough everyday lab work, you need a reliable and precise magnetic plate system to coordinate your articulators. Adesso Split magnetic plate systems and matching calibration keys are available for practically all well-known articulators. Adesso Split not only enables articulators with the same construction to be calibrated with each other, but also allows certain different articulators to be coordinated with one another, as the illustration to the right shows. Do you also own articulators from various manufacturers and are unsure whether it is possible to coordinate them? Ask us - we will be happy to help you further!



...do you still have questions?
We look forward to your call...
Fon +49 7236-9 33 69 0

Simple troubleshooting

Is calibration still intact? The **adesso split® calibration key** provides you with the answer.



Using a divided **Adesso Split calibration key**, a defect in the articulator or a fault in coordination can be identified immediately and with ease

The **Adesso Split calibration key** is an adjustment and testing tool in one. Thanks to the patented construction, the key can be divided in two for testing purposes. In contrast to one-piece calibration keys, the **Adesso Split calibration keys** enable even the smallest inaccuracies to be identified quickly and easily. Even defective condyle housings or distorted frames on your articulators can be identified quickly using the **Adesso Split calibration key**.



adesso split® calibration key

The original Adesso Split calibration keys from Baumann-Dental are available for all common articulator symmetries. The product is supplied in a practical box. Are you unable to find what you are looking for here, or do you still have questions? Give us a call!

10025	Adesso Split calibration key Artex ¹ , Arto 116 mm, Stratos, Reference
10025-C	Adesso Split calibration key Artex ¹ Carbon, Arto 126 mm, Physio Logic
10027	Adesso Split calibration key / Denar
10030	Adesso Split calibration key / KaVo ³
10031	Adesso Split calibration key / Panadent 108 mm
10031-01	Adesso Split calibration key / Panadent 93 mm
10032	Adesso Split calibration key / Rational
10033	Adesso Split calibration key / SAM ² P
99999-01	Box for calibration key / Adesso Split
10024	Allen key / 5 mm / for calibration key



Tip: Naturally, you can integrate Adesso Split yourself and calibrate your articulators on site. If you would prefer to leave this to our workshop experts, that's no problem, and we'll be happy to take care of it for you. You will be surprised how affordable our workshop's work is.




Figure shows: Parts supplied with an **Adesso Split calibration key**



CONSISTENT, CROSS-SYSTEM CALIBRATION

The patented Adesso Multisplit magnetic plate system perfectly combines the two most common magnetic plate systems with no loss of precision: the proven Adesso Split system from Baumann-Dental and the Splitex¹ system from AmannGirrbach.



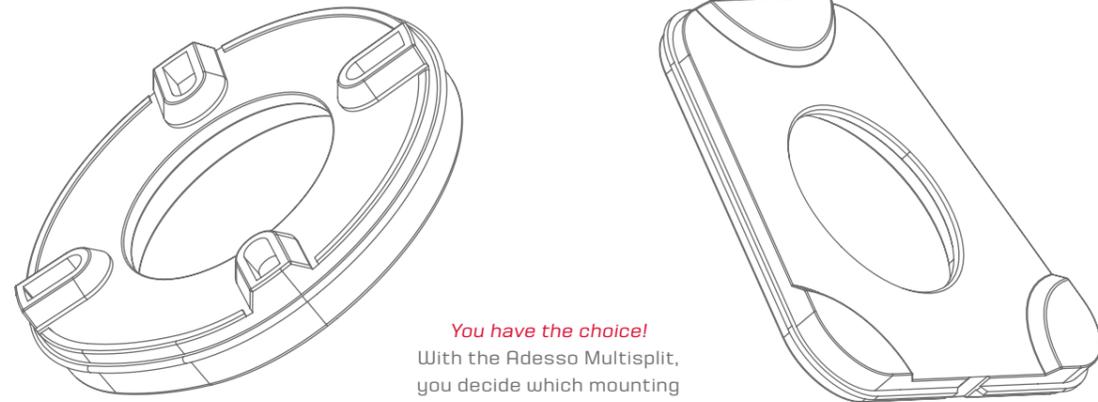

You have the choice!

With the Adesso Multisplit, you decide which mounting plates you work with: whether the original Adesso Split mounting plates from Baumann-Dental or the mounting plates for the Splitex¹ system from AmannGirrbach.



Adesso MULTISPLIT®

The multifunctional magnetic plate system

With the patented Adesso Multisplit magnetic plate system from Baumann-Dental, you can use both of the most popular magnetic plate systems in one articulator, or scanner, or mounting articulator, or transfer stand, etc., etc.



Your benefits:

- Cross-system calibration of articulators
- Mounting plates for the systems Adesso Split and Splitex¹ can optionally be used
- Other accessories from Adesso Split and Splitex¹ have unlimited compatibility
- Adesso Multisplit can be calibrated with the centring keys for the Adesso Split or Splitex¹ systems
- Low costs through dual benefits

Figure shows Artex¹ Carbon articulator with Adesso Multisplit set #10200-C and spacer plate #10265, for reducing the Carbon articulator to the construction height of 116 mm.

*...do you still have questions?
We look forward to your call...
Fon +49 7236-9 33 69 0*





Figure shows **Adesso Multisplit transfer plate #10200-15** for 3Shape® scanner. The extremely flat transfer plate makes it possible for the first time to directly place models on Adesso Split and Splitex® bases directly in the scanner. Not suitable for 3Shape® D900L scanners!



adesso MULTISPLIT® complete systems

Multifunctional magnetic plate system for calibration.

10200	Adesso Multisplit Set / Artek®, Arto
10200-C	Adesso Multisplit set / Artek® Carbon / without magnets

adesso MULTISPLIT® individual components *without magnets*

10200-10	Adesso Multisplit base plate / Artek®, Arto
10200-11	Adesso Multisplit adjusting base plate / Artek®, Arto
10200-15	Adesso Multisplit transfer plate for 3Shape® scanner
10200-C-10	Adesso Multisplit base plate / Artek® Carbon
10200-C-11	Adesso Multisplit adjusting base plate / Artek® Carbon
10200-C-12	Adhesive plate Adesso Multisplit / Splitex® for set Artek® Carbon
10265	Adesso Multisplit spacer plate 10 mm e.g. for Artek® Carbon without illustration
10266	Adesso Multisplit spacer plate 20 mm for Physio-Logic without illustration
10020	Adesso Split magnet / without illustration / graduated price! for use with Adesso Split, Adesso Multisplit plates and other magnetic plate systems, e.g. Splitex® or Quicksplit®



Figure shows: Parts supplied with **Adesso Multisplit Artek Carbon #10200-C**

The newly designed milling base in the version adesso MULTISPLIT®

Particular value has been placed on an extremely flat construction in order to process higher models including base. You can determine the working height yourself beforehand. The milling base is exceptionally stable and easy to use.

The features in detail:

- Extremely flat construction
- Can be rotated in several directions by up to 90°
- Adjustable height
- Magnetic base plate
- Also available for Adesso Split and Model System 2000

10039-02	Adesso Multisplit milling base
----------	--------------------------------



Figure shows: **Milling base** with Adesso Multisplit base plate

Splitex¹ compatible magnetic plate system

The tried-and-tested plate system at an especially attractive price. Available for Artex¹, Arto and Artex¹ Carbon.

Splitex¹ compatible complete systems

10010	Splitex ¹ compatible magnetic plate set / Artex ¹ / Arto
10010-C	Splitex ¹ compatible magnetic plate set / Artex ¹ Carbon / without magnets

Splitex¹ compatible single components *without magnets*

10010-10	Splitex ¹ compatible base plate / Artex ¹ / Arto
10010-11	Splitex ¹ compatible adjusting base plate / Artex ¹ / Arto
10010-C-10	Splitex ¹ compatible base plate / Artex ¹ Carbon
10010-C-11	Splitex ¹ compatible adjusting base plate / Artex ¹ Carbon
10066	Splitex ¹ compatible spacer plate 10 mm for Artex ¹ Carbon / with magnet
10020	Adesso Split magnet / without illustration / graduated prices! for use with Adesso Split, Adesso Multisplit plates and other magnetic plate systems, e.g. Splitex ¹ or Quicksplit [®]



Figure shows parts supplied with the Splitex¹ compatible plate set for Artex¹ #10010

Figure shows parts supplied with the Splitex¹ compatible plate set for Artex¹ Carbon #10010-C



Figure shows: Artex¹ Carbon from AmannGirrbach with Splitex¹ compatible plate set for Carbon #10010-C

...do you still have questions?
We look forward to your call...
Fon +49 7236-9 33 69 0

Figure shows parts supplied with magnetic plate for the system **KaVo³-Split** #99997-10. The plate can be mounted with the enclosed mounting material in the upper part of the articulator. Alternatively, it can be adhered to the lower part.



Tip: In order to achieve highly precise and, above all, reproducible calibration in your valuable KaVo³ articulators, we recommend integrating our new base plate in the upper part of the articulator. To reduce costs, the available original KaVo³ plate can be fitted in the lower part of the articulator. The plates are also suitable for integration in scan articulators and fixators. Give us a call...



KaVo³-Split compatible magnet plate & centring key

People who own KaVo³ articulators used to be instructed to send their articulators to the factory to check the coordination. This was time-consuming and expensive.

The check was hampered by the fact that both magnetic plates in the KaVo³-Split system are suspended 'floating'. We took a close look at this issue. Using our proven Adesso Split calibration key as a role model, we have designed a divisible calibration key for the KaVo³-Split system, which provides you with the perfect adjustment and testing tool to verify your valuable KaVo³ articulators in your lab. The specially-made base plate for the KaVo³-Split system is also new.

99997-10	KaVo ³ -Split compatible base plate / incl. magnet
99998	Calibration key for KaVo ³ -Split system
10020	Adesso Split magnet / without illustration / graduated prices! for use with Adesso Split, Adesso Multisplit plates and other magnetic plate systems, e.g. Splitex [®] or Quicksplit [®]
99999-01	Box for calibration key / KaVo ³ Split

*...do you still have questions?
We look forward to your call...
Fon +49 7236-9 33 69 0*



Model mounting

Plastic mounting plates for the following systems: Adesso Split®/ Splitex¹ / KaVo³ Split

The right one for everyone...

For the perfect mounting of your models, Baumann-Dental offers you a variety of different mounting plates with which we serve the most common magnetic plate systems on the market. With the outstanding quality you expect. Naturally, all of the products come from our own production facility. By building new tools and through improved production methods, we have been able to significantly boost the precision of the mounting plates once again and you will be delighted to know that our prices are falling!

Please also pay particular attention to our fabulous new Connect mounting plates, which make it practically impossible for the plates to inadvertently loosen from the plaster base.



- Perfect hold in the plaster
- Easy to clean
- Can be reused multiple times
- Extremely flat construction
- Integrated magnetic counter



Figure shows: *Connect mounting plates* for the systems *Adesso Split* and *Splitex*¹

CoNECT mounting plates adesso split®

14050	Connect for Adesso Split / grey / 1 pc.
1405020	Connect for Adesso Split / grey / 20 pcs.
1405050	Connect for Adesso Split / grey / 50 pcs.
14051	Connect for Adesso Split / white / 1 pc.
1405120	Connect for Adesso Split / white / 20 pcs.
1405150	Connect for Adesso Split / white / 50 pcs.

CoNECT silicone matrix adesso split®

11004	Silicone matrix for Connect
-------	-----------------------------

CoNECT mounting plates Splitex¹

14030	Connect for Splitex ¹ / white / 1 pc.
1403020	Connect for Splitex ¹ / white / 20 pcs.
1403050	Connect for Splitex ¹ / white / 50 pcs.



Image shows *Connect mounting plate* from the underside. *Connect* does not require any adhesive disk as this is already firmly integrated.

! *Tip: For particularly clean and attractive mounting bases, suitable matrices are available for the Adesso Split Connect plates. These are simply placed around the mounting plate before insertion in the articulator and removed again when the base is removed.*

CoNECT® The brilliant plaster claw!
for the systems *adesso split®* and *Splitex*¹

In modern dental technology, low-expanding articulation plasters are increasingly being used, which certainly benefits precision. As a result of the low expansion values, mounting plates inadvertently loosen out of the plaster base more and more frequently. In addition, more and more high model systems are being used, yet space in the articulator is limited. So, what should be done? This discussion has been brought to an end with the market launch of the new Connect mounting plates: The plates are very flat,

thus providing you with additional space. Through the innovative retention surface, there are more undercuts than necessary, and your models are fixed completely secure in the base. The innovative Connect plates bring with them the adhesive disks - which is also very practical - and, not only that, but they are also easy to clean and can be used multiple times. Request your no-obligation test sample now!

Connect mounting plates from Baumann-Dental - always a step ahead!

CoNECT Insert mounting plate in the articulator



CoNECT Apply matrix and fill with plaster



Place and orient the model on the *CoNECT* plate



**YOU'RE DONE!
IT REALLY IS THAT SIMPLE!**

'Basic' plastic mounting plates

Using our plastic mounting plates guarantees optimum precision and cost-effectiveness in the lab. The plates enable clean model mounting and can be reused numerous times. The product is supplied without adhesive disks.



Mounting plates for Adesso Split® Choose between 'basic' and 'snap-in' versions

In the basic version, the adhesive disk is simply put in place and then held in position by the articulation plaster. With the snap-in version, the adhesive disk is secured in the mounting plate before plastering using small claws.

Plastic mounting plate / units per pack 100 pcs. / delivery without adhesive disk / matching adhesive disks: #10028 or #11010 (page 59)

14100	Adesso Split mounting plate / yellow / basic
14100.2	Adesso Split mounting plate / yellow / snap-in
14200	Adesso Split mounting plate / blue / basic
14200.2	Adesso Split mounting plate / blue / snap-in
14300	Adesso Split mounting plate / white / basic
14300.2	Adesso Split mounting plate / white / snap-in
14400	Adesso Split mounting plate / mint green / basic
14400.2	Adesso Split mounting plate / mint green / snap-in
14500	Adesso Split mounting plate / red / basic
14500.2	Adesso Split mounting plate / red / snap-in
14900	Adesso Split mounting plate / orange / basic
14920	Adesso Split mounting plate / grey / basic
14931	Adesso Split mounting plate / black / basic

Mounting plates system SpliteX¹

Plastic mounting plate / units per pack 100 pcs. / delivery without adhesive disk / matching adhesive disks: #10028 or #11010 (page 59)

14002	Mounting plate / basic / for SpliteX ¹ / yellow
14003	Mounting plate / basic / for SpliteX ¹ / blue
14941	Mounting plate / basic / for SpliteX ¹ / black
14942	Mounting plate / basic / for SpliteX ¹ / white



Mounting plates system KAVO³- Split

Plastic mounting plate / 50 pcs. / delivery without adhesive disk / matching adhesive disk #14028 (page 59)

14020	Mounting plate / basic / for KaVo ³ Split / white
14020-10	Mounting plate / basic / for KaVo ³ Split / blue
14020-11	Mounting plate / basic / for KaVo ³ Split / red
14020-12	Mounting plate / basic / for KaVo ³ Split / green



You can find suitable adhesive disks on the following page >>>



Adhesive disks
for plastic mounting plates *diameter 35 mm*

e.g. for type Adesso Split / Splitex¹ / Quicksplit⁵ / KaVo³ /
Please ask us for our graduated prices.

10028	Adhesive disks / galvanic finish / 100 pcs.
11010	Adhesive disks / magnetic stainless steel / 100 pcs.
14028	Adhesive disks for KaVo ³ mounting plates / 50 pcs.



The adhesive disks for **KaVo³** mounting plates have a layer so they can be applied flush to the mounting plate. This increases adhesion strength enormously.



Mounting plates for adesso split[®]
Premium version - the special one!

Adhesive strength can be infinitely regulated with the premium mounting plate.

14006	Mounting plate for Adesso Split Premium / white / 50 pcs. / delivery incl. counter screws
-------	---

Stock remainders are available in other colours. Please ask us!



The simple Allen key #14005 can be used to adjust the integrated counter screw; the adhesion strength can be easily adjusted to suit your needs.

...do you still have questions?
We look forward to your call...
Fon +49 7236-9 33 69 0

Model Making

The right system for **everyone!**

No matter whether you prefer classic plaster models, models on premade model plates according to the principle of the Zeiser⁴ system, Combiflex or our proven pin base plates from the Model System 2000 or Quadra Base.

At Baumann-Dental, you too can find the right system for your needs. Quality and cost-effectiveness take priority for all model systems from Baumann-Dental.

Let yourself be inspired...

 **Modellsystem 2000**®
Baumann-Dental

 **Duett System**®
Baumann-Dental

 **comBi flex**®
Baumann-Dental

 **Quadra Base**®
Baumann-Dental



**Cost-effective:***Can be reused up to 50 times**70% time saving**No machines required***Precise:***No multiple plaster expansions**Segment sits without wobbling**Accurately fitting splitcast***Ergonomic:***No drills or tweezers required**Easy sawing without damaging pins**Optimum checking of fit***Model System 2000**®

The gold standard when it comes to saw cut models..

In just a few minutes, the Model System 2000 can produce a high-quality saw cut model, which offers outstanding precision that results from just one expansion process.

The optimised design and the sophisticated arrangement of the pins guarantees ergonomic setting of the individual segments. Even critical approximal contacts can be precisely set on the model.

The already integrated splitcast is characterised by an extremely flat construction in order to facilitate model mounting in the articulator. Individual pins can be easily removed where required in order to create space for a model analogue, for instance.



Model System 2000® starter kit

Pin base plate / large / 2 pcs. / #19000

Pin base plate / small / 3 pcs. / #19020

Matrix / large / 1 pc. / #19060

Matrix / small / 2 pcs. / #19061

Archiving plate Eco / large / 2 pcs. / #19050

Archiving plate Eco / small / 3 pcs. / #19051

Positioning plate / large / transparent / 1 pc. / #19036-01

Positioning plate / small / transparent / 1 pc. / #19035-01

Model segment removal key / large / 1 pc. / #19030

Table for model segment removal key / small / 1 pc. / #19033

Separating fluid / spray bottle / 100 ml / 1 pc. / #19080

19090 Starter kit / composition see above

19091 Starter kit incl. cleaning spray / #19088

19090/E Starter kit incl. Quadra Base

Figure shows: Parts supplied
 starter kit MS 2000 #19090

Model System 2000® Components

Pin base and splitcast plates

Graduated prices from 10 pcs.

19000	Pin base plate / large
19020	Pin base plate / small
19100	Pin base plate / large / without splitcast plate
19111	Pin base plate / small / without splitcast plate
19110	Splitcast plate / large
19112	Splitcast plate / small

Matrices

19060	Matrix / large
19061	Matrix / small

Archiving plates

19050	Archiving plates Eco / large / 25 pcs.
19051	Archiving plates Eco / small / 25 pcs.
19050-DE LUXE	Archiving plates / deluxe / large / 25 pcs.
19051-DE LUXE	Archiving plates / deluxe / small / 25 pcs.

Insulating and cleaning

19080	Separating fluid / spray bottle / 100 ml
19081	Separating fluid / refill canister / 500 ml
19088	Cleaning spray / 400 ml
19089	Reamer

Other accessories

19030	Model segment removal key / large
19031	Model segment removal key / small
19032	Table for model segment removal key / large
19033	Table for model segment removal key / small
19001	Model removal plate / large
19021	Model removal plate / small
19036-01	Positioning plate large / transparent / 2 pcs.
19035-01	Positioning plate small / transparent / 2 pcs.





Model System 2000® special cuvette

Special cuvette for the Model System 2000 for the duplication of individual or numerous tooth segments in the correct position for individual ceramic layering.

Thanks to the inner divider, which can be used flexibly, segments can be duplicated in a way that saves material.

In the future, work with investment material segments directly on your Model System 2000®, without laborious translation to a conventional investment material system.

Content:

- 1 cuvette for pin base plate small
- 3 inner dividers
- 1 model segment removal key
- 1 instruction manual

62210 MS 2000 cuvette / small

*...do you still have questions?
 We look forward to your call...*
Fon +49 7236-9 33 69 0

Model System 2000® milling base

The new milling base for the Model System 2000. Particular value has been placed on an extremely flat construction in order to process higher models including base. You can determine the working height yourself beforehand. The milling base is exceptionally stable and easy to use. The base plate can hold large and small pin base plates, which are fixed completely secure using strong magnets.

The features in detail:

- Extremely flat construction
- Can be rotated in several directions by up to 90°
- Adjustable height
- Magnetic base plate
- Also available for Adesso Split and Adesso Multisplit

10039-03 Milling base Model System 2000



Duet plate system

The Duet plates are used in all areas of dental technology and dentistry for the capping and mounting of all unsawn models. The plates can be reused numerous times and a splitcast plate is integrated.

Application e.g. for:

- Counterbites
- Situation models
- Models for functional diagnostics and treatment
- Splint therapy
- Models for combination and milling work
- Orthodontic models (set-up)
- Models for total prosthetics
- Backing
- Implantology models, etc.



Using a *different colour code*, the various models and situations can not only be linked together in the functional diagnoses, but all models can also be exchanged for one another. Both parts of the Duet plate are thus produced in different colours, each with a different meaning: (e.g. in functional diagnostics). If treatment is provided via rearticulation, the colour of the splitcast plate indicates the stage, while the model remains the same. If treatment is provided via diagnostic models, the colour of the model base plate indicates the stage, while the splitcast plate remains the same. The Duet plates can be reused.



*...do you still have questions?
 We look forward to your call...*
Fon +49 7236-9 33 69 0



Initial situation for analysis and planning.



Therapeutic position here with splints.



Centric position with diagnostic wax-up.



Lower jaw situation prepared in def. centric against upper jaw situ.



Lower jaw situation prepared in def. Centric against upper jaw wax-up.

Duet plate system

19095	Duet plates starter kit / complete
19070	Duet plate / large / white
19071	Duet plate / large / red
19072	Duet plate / large / blue
19074	Duet plate / large / green
19060	Matrix / large

Treatment set Parts supplied

White plates / large / 4 pcs. / #19070
Blue plates / large / 4 pcs. / #19072
Red plates / large / 4 pcs. / #19071
Green plates / large / 4 pcs. / #19074
Sleeves / large / 2 pcs. / #19060



Please note: The 'small' Duet plates are no longer manufactured. Some stock remainders are still available; please enquire if interested.



Duet System[®]
Baumann-Dental





The plastic version of the proven Model System 2000 with the precision you expect from Baumann-Dental at an attractive price. The Evo system impresses with its flexibility and simplicity.

The new Model System 2000 Evo pin base plate is made completely from high-quality plastic.

The Evo pin base plate corresponds to the size of our classic MS 2000 pin base plate small #19020.



Model System 2000® EVO Introduction set

Pin base plate / small / 10 pcs. / #1970010

Insert with metal core / 10 pcs. / #1961010

Matrix / small / 2 pcs. / #19061

Splitcast plate / small / 2 pcs. / #19112

Separating fluid / spray bottle / 100 ml / 1 pc. / #19080

Hammer for model systems - plastic head / 1 pc. / #11102-02

19690 Model System 2000 Evo starter kit

Model System 2000® EVO components

1970010	Pin base plates / small / 10 pcs.
1961010	Insert with metal core / 10 pcs. /
19620/N10	Insert without metal core / 10 pcs.
19112	Splitcast plate small / 1 pc. / white
19080	Separating fluid / spray bottle / 100 ml
19081	Separating fluid / refill canister / 500 ml
11102-02	Hammer for model systems - plastic head
19061	Matrix / small / 1 pc.



MS 2000 Evo Flexibility even for model mounting...



Evo pin base plates from Baumann-Dental are constructed in such a way that they can also be articulated directly, without a splitcast plate, making them an extremely cost-effective version.

For this you need:

- #1970010 pin base plate / small
- #19620/N10 insert for pin base plate Evo

The extremely precise yet very cost-effective **MS 2000 Evo** pin base plates are specially suited to large dental labs, as they were developed for one-time use.

Also suitable for printing and milling techniques, the **MS 2000 Evo** pin base plate is now supported by software applications from major manufacturers.



The **MS 2000 Evo system** is also an affordable alternative when you prefer the classic method of the 'double splitcast' when mounting models. You need:

- #1970010 pin base plate - small
- #1961010 insert for pin base plate with metal core
- #19112 Splitcast plate - small



comBi flex®

Base plate premium with pressed-in magnetic counter

The Combiflex model system from Baumann-Dental is a precise and affordable alternative. We offer you a complete system with all the necessary individual components, suitable for the pin drilling units from Zeiser¹ and AmannGirrbach. The short basal positioning elements in the Combiflex plates prevent lateral sliding of the model base on the splitcast plate.

Are you interested and would like to try Combiflex with no obligation? No problem, just get in touch with us!

Combiflex, the alternative...



*proven millions of times over
outstanding price-performance ratio
pressed-in magnetic counter
top precision **made in Germany***

Premium base plate **large** / units per pack 100

With pressed-in magnetic counter. (Corresponds to Giroform¹ size XL)

11115	Base plates premium / large / white
11115-10	Base plates premium / large / blue



Premium base plate **small** / units per pack 100

With pressed-in magnetic counter. (Corresponds to Giroform¹ size L)

11116	Base plates premium / small / white
11116-10	Base plates premium / small / blue
11116-11	Base plates premium / small / terracotta
11116-12	Base plates premium / small / black




*...do you still have questions?
We look forward to your call...
Fon +49 7236-9 33 69 0*



comBi flex® base plate Basic / units per pack 100 pcs.

Variant of the Combiflex system. For use with counter screw M8.2 #10059 or with the Combiflex metal pins #10085 ('large' = Giroform¹ XL, 'small' = L)

11109	Base plates basic / large / white
11110	Base plates basic / small / white
11110-10	Base plates basic / small / blue

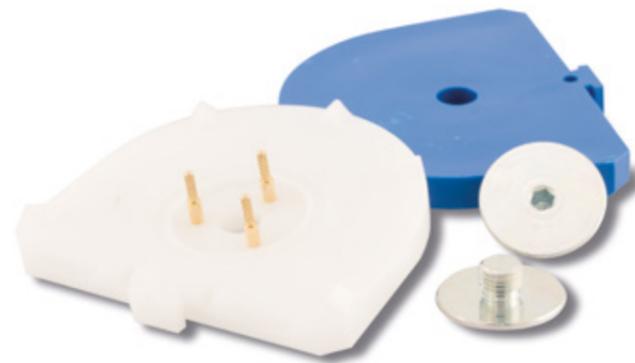

comBi flex® splitcast plate

The fast way to achieve a perfect splitcast. Secures the base plate perfectly against lateral displacement. The sophisticated construction supports the base plates in all areas.


1111825PK	Splitcast plate - large, 25 pcs.
1111850PK	Splitcast plate - large, 50 pcs.
1111925PK	Splitcast plate - small, 25 pcs.
1111950PK	Splitcast plate - small, 50 pcs.



Plate holders

Our workshop service modifies available plate holders e.g. from Zeiser⁴ or Giroform¹ to hold our Combiflex plates. This service is free of charge for those switching from other systems.

176000	Combiflex universal plate holder for Zeiser ⁴ P10 and Giroform ¹ pin drilling units / fits large and small Combiflex, Zeiser ⁴ and Giroform ¹ base plates.
176000-01	Locking screw for universal plate holder
176112	Plate holder for Zeiser ⁴ / small / black / without illustration
176122	Plate holder for Zeiser ⁴ / large / blue
176142	Combiflex plate holder, large for Zeiser ⁴ and AmannGirrbach pin drilling units of older construction types.

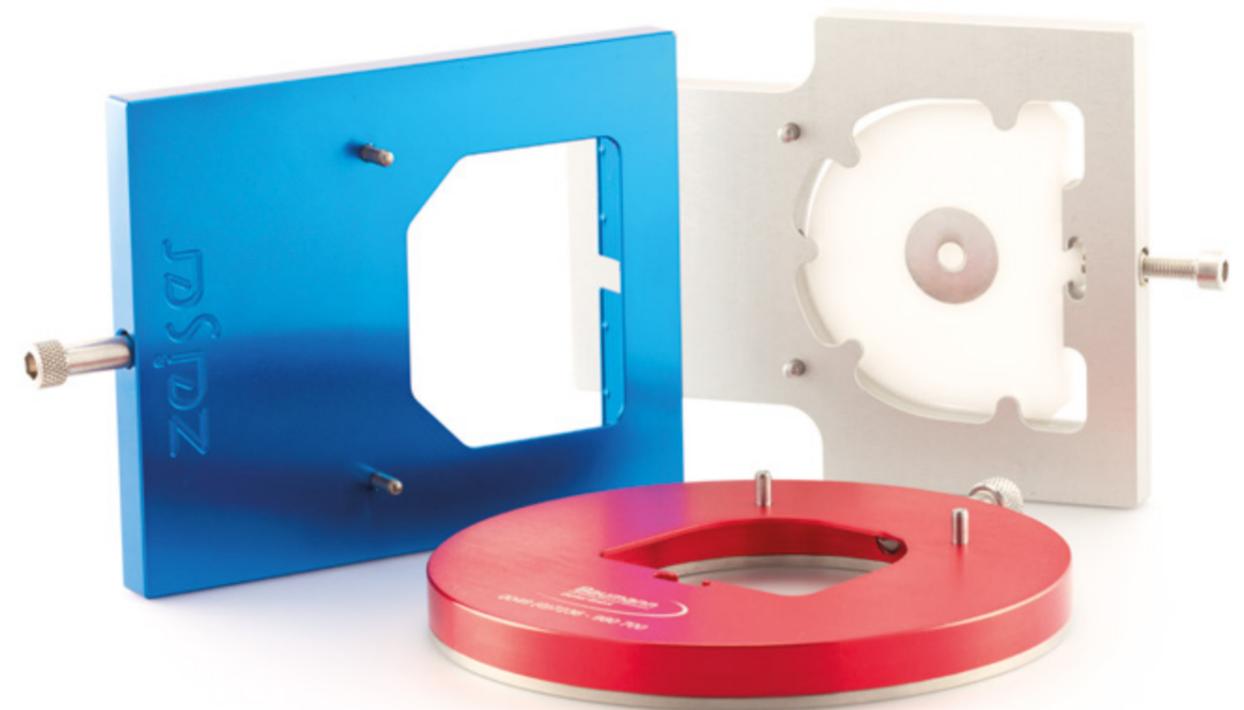

COMBI FLEX® accessories for Zeiser⁴ and Giroform¹

We offer you the full range of accessories for the systems Combiflex, Zeiser⁴, Giroform¹ and similar systems.

Are you unable to find what you are looking for here?

Give us a call, we'll be happy to help you further...

10085	Pins 14 mm for Combiflex / Zeiser ⁴ / Giroform ¹ and similar / 1,000 pcs. / Please ask us for our graduated prices.
10085-k	Pins compact 12.5 mm for Combiflex / Zeiser ⁴ / Giroform ¹ and similar / 1,000 pcs. / Please ask us for our graduated prices.
10082	Pin drill conical for Combiflex / Zeiser ⁴ / Giroform ¹ and similar
10059	M8.2 screw for Combiflex, Zeiser ⁴ , Giroform ¹ / 100 pcs. / for use with our 'basic' base plates
14005	T-handle Allen key for use with our M 8.2 screws
11107	Combiflex putty / pliable silicone / silver-grey / 1 kg / for Combiflex / Zeiser ⁴ / Giroform ¹ and similar
11170	comBiflex Putty natural / natural pliable material / 1kg for Combiflex / Zeiser ⁴ / Giroform ¹ and similar



Tip: The compact 12.5 mm pin has the same shaft as the large 14 mm pin, so there is no need to change drill bit. This pin is very well-suited to the manufacture of very small or narrow lower jaw models, as it has a shorter and narrower retention surface.



Positioning plates for COMBI FLEX® / Giroform¹ and similar

Facilitates the positioning of impressions on the impression carrier


11111	Positioning plate / transparent / large / 2 pcs.
11112	Positioning plate / transparent / small / 2 pcs.



Impression carrier for COMBI FLEX® and Giroform¹

The impression carriers can be used for all base plates and fit the plate holders and drilling units from Zeiser⁴ and AmannGirrbach. Impression carrier L is 75 mm high, making it very well-suited to holding even high, heavily offset impression trays and positioning these in the correct place.

11120	Combiflex impression carrier standard / size M / 50 mm
11121	Combiflex impression carrier / size L / 75 mm

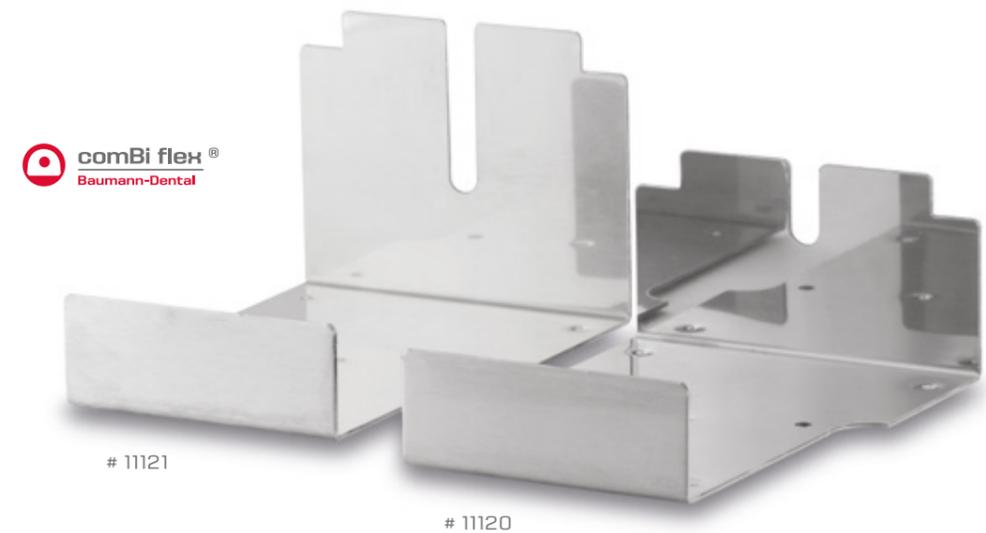



...do you still have questions?
We look forward to your call...
Fon +49 7236-9 33 69 0

Zeiser⁴ base plates

Top quality at a low price

Available in the versions 'premium', with pressed-in magnetic counter, and 'basic', for use with our magnetic counter screw #10059 or with pins #10085 / #10085-K

Premium / units per pack 100 pcs.

11105	Base plate for Zeiser ⁴ premium / large / var. B / white
11106	Base plate for Zeiser ⁴ premium / small / var. B / white

Basic / units per pack 100 pcs.

11103	Base plate for Zeiser ⁴ basic / large / white
11103-10	Base plate for Zeiser ⁴ basic / large / blue
11104	Base plate for Zeiser ⁴ Basic / small / white
11104-10	Base plate for Zeiser ⁴ Basic / small / white
10059	M8.2 screw for Combiflex, Zeiser ⁴ , Giroform ¹ / 100 pcs.
14005	T-handle Allen key for M8.2 screw / 1 pc.

Positioning plate

11113	Positioning plate for Zeiser ⁴ / transparent / large / 2 pcs. / without illustration
11114	Positioning plate for Zeiser ⁴ / transparent / small / 2 pcs. / without illustration

Other accessories for model manufacture according to Zeiser⁴, such as the pins shown can be found in the 'Combiflex' section

Our Zeiser⁴ type base plates offer *impressive quality* at a low price. They are not suitable for use with the original Zeiser⁴ splitcast plates.





Quadra Base in use with the practical disposable articulator.

Quadra Base® Model manufacture reduced to the essentials...

The Quadra Base model system from Baumann Dental is a pure quadrant model system. It is based on the principle of our proven Model System 2000. Quadra Base plates can be articulated or used with our practical disposable articulators. When using quadrant impressions, Quadra Base is an affordable and fast alternative. Try it for yourself...

Quadra Base® introduction set

- Model base right / 10 pcs. / #19200-1
- Model base left / 10 pcs. / #19200-2
- Matrix right / 2 pcs. / #19260-1
- Matrix left / 2 pcs. / #19260-2
- Disposable articulator / 10 pcs. / #19270
- Model adapter / 20 pcs. / #19271
- Separating fluid / spray bottle / 100 ml / #19080

19290	Quadra Base starter kit
-------	-------------------------

single components

19200-1	Model base right / 50 pcs.
19200-2	Model base left / 50 pcs.
19270	Disposable articulator / 10 pcs.
19271	Model adapter / 20 pcs.
19260-1	Matrix right / 1 pc.
19260-2	Matrix left / 1 pc.
19080	Separating fluid / spray bottle / 100 ml



Cost-effective
Time saved when manufacturing models
No machines required

Precise
No multiple plaster expansions
Segment sits without wobbling

Ergonomic
No drills or tweezers required
Easy sawing without damaging pins
Optimum checking of fit



Model base moulds

Firmly integrated magnet and supplied without matrix.
Some of the matching matrixes are available in different heights - please select the one you require.

small

10072	Model base mould / small / grey
10076	Matrix for model base mould / small / 21.7 mm
10073	Matrix for model base mould / small / 23.5 mm
10077	Matrix for model base mould / small / 27.0 mm

large

10070	Model base mould / large / blue
10078	Matrix for model base mould / large / 21.7 mm
10071	Matrix for model base mould / large / 23.5 mm
10079	Matrix for model base mould / large / 27.0 mm

extra large

10074	Model base mould / extra large / black
10075	Matrix for model base mould / extra large / 23.5 mm



Figure shows the cast-in magnets in the model base mould large #10070



! Tip: Baumann-Dental model base moulds also fit our Xpress pin system perfectly!

Xpress pins

For lovers of the classic plaster model

The Xpress pin system from Baumann-Dental is extremely flexible to use. The pins can be fixed in place in the tothing without adhesive. The special design means the pins slide cleanly into the sleeves and ensure a perfect hold. Optionally available approximal pins secure the segments against mesio-distal pull-out torque and optimise the plaster model for the best possible setting of the approximal contacts. Ranges with closed and open sleeves are available.

10810	Brass pins incl. plastic sleeves / closed / 1,000 pcs.
10813	Brass pins incl. plastic sleeves / open / 1,000 pcs.
10811	Approximal pins / 1,000 pcs.
10812	Xpress pin system / drill / 1 pc. / without illustration



The Xpress pins are inserted into the tothing by the conical side and are fixed in place by tapping gently with a small hammer. There is no need to use adhesive on the pins.

The optionally available approximal pins are inserted into the holes in the tothing in such a way that the milled recess remains visible. The recess is later anchored in the base plaster.

The Xpress pin sleeves are available open or closed. These give the pins a secure, wobble-free hold.



Figure shows: Xpress pins and approximal pins.

Model manufacture Splitcast item / accessories

Accessories

16250	Occlu-contact spray / 75 ml / blue
16251	Occlu-contact spray / 75 ml / green
11009	Silicone care and insulation spray / 400 ml
20000-10PK	Membrane box / 36x36x16 mm / 10 pcs.

Our **silicone spray** #11009 is perfectly suited to caring for your valuable articulators, magnetic plate systems, etc. It is also suitable for insulating your Adesso Split master plates.



16250 / 16251



20000-10PK



Splitcast items

10052	Retentive discs / punched / galvanic finish / 100 pcs.
10053	Magnet housings pressed / galvanic finish / 50 pcs.
10054	Retentive discs / milled / galvanic finish / 100 pcs.
10055	Magnet housings milled / extra flat / 25x5 mm galvanic finish / 50 pcs.
10056	Power magnets 14x3 mm / 100 pcs.
10055/56	Magnet housings milled / extra flat incl. power magnet / set / 10 pcs.
10057	Magnets round 20x6 mm / 100 pcs.
10058	Magnet housings milled / galvanic finish / 50 pcs.
10028	Adhering discs / galvanic finish / 100 pcs.
11010	Adhering discs / magnetic stainless steel / 100 pcs.



...do you still have questions?
We look forward to your call...
Fon +49 7236-9 33 69 0

Grinding Equipment

Abrasive bands for dry trimmer

The grinding belts from Baumann-Dental are available for the most common dry trimmers and in various grit sizes. The optimised surface ensures high abrasion with a very long service life. Sales unit: 10 pcs.

Harnisch & Rieth / 100 x 880 mm

25000-40	Grinding belt for Harnisch & Rieth / grit 40
25000-50	Grinding belt for Harnisch & Rieth / grit 50
25000-60	Grinding belt for Harnisch & Rieth / grit 60
25000-120	Grinding belt for Harnisch & Rieth / grit 120

KaVo³ / 100 x 950 mm

25010-40	Grinding belt for KaVo ³ / grit 40
25010-60	Grinding belt for KaVo ³ / grit 60
25010-80	Grinding belt for KaVo ³ / grit 80

AmannGirrbach / 90 x 610 mm

25020-60	Grinding belt for AmannGirrbach / grit 60
25020-100	Grinding belt for AmannGirrbach / grit 100



Grinding belt mandrels & grinding belts

can be used in the technician hand piece, shaft: 2.35 mm

Ringband shape / 13 x 12 mm

25300	Mandrel 2,35 mm for grinding belts ringband shape / 1 pc.
25300-80	Grinding belts ringband shape 13x12 mm / grit 80 / 50 pcs.
25300-120	Grinding belts ringband shape 13x12 mm / grit 120 / 50 pcs.

Bullet shape

25400	Mandrel 2,35 mm for grinding belts bullet shape / 1 pc.
25400-80	Grinding belts bullet shape / grit 80 / 10 pcs.

Tungsten Carbide cutter / 6 mm

25202	Tungsten carbide cutter conical / 6 mm shaft / e.g. for AmannGirrbach Multitrim
-------	---

...do you still have questions?
We look forward to your call...
Fon +49 7236-9 33 69 0



Figure shows: *Duplicator* #62100 with spacer bushing. The Adesso Multisplit plate #10200-10 shown is not supplied.

Duplicating technique

More modern than ever before -the verticator principle

The duplicator from Baumann-Dental was developed as a universal and highly precise vertication tool. The field of application is practically unlimited: whether for transferring a wax-up to the master model, for general gnathological work to transfer the occlusion, the splint technique, for non-prep partial crowns or even for linings and remountings. The duplicator is supplied with different spacer bushings to adapt the tool individually to the respective work requirements. The standardised mountings in the lower part also enable magnetic plates to be easily installed.

62100	Duplicator complete
62100-01	Duplicator / upper part

Basic duplicating cuvette

Duplicating system for highly precise and cost-effective duplications

- Simple and quick to use
- The duplicating material adheres excellently in the former
- Extremely cost-effective
- Very well-suited to long-term archiving, patient data can be directly noted on the former
- Stable, long-lasting cuvette base made from anodised aluminium

61000	Basic duplicating cuvette / starter kit / cuvette base & 12 formers
61100-01	Basic duplicating cuvette / former / 1 pc.
61200	Basic duplicating cuvette / cuvette base / 1 pc.



...do you still have questions?
 We look forward to your call...
Fon +49 7236-9 33 69 0

Transfer equipment

Transfer stands & transfer jig

Transfer stand universal

For transferring the Artex¹ face bow registration in the correct position in an Arto, Artex¹, Carbon or SAM² articulator (in conjunction with the adapter plate for SAM² #10046). For retrofitting with a base plate that can be screwed in place or for use without a base plate, e.g. for a transfer table that can be screwed in place.

10049 Transfer stand universal

adesso Split[®] Transfer stand

For transferring the Artex¹ face bow registration in the correct position in an Arto, Artex¹ or Carbon articulator, incl. Adesso Split base plate, for use with our Adesso Split transfer table (#10048).

10050 Adesso Split transfer stand

adesso MULTISPLIT[®] Transfer stand

For transferring the Artex¹ face bow registration in the correct position in an Arto, Artex¹ or Carbon articulator, incl. Adesso Multisplit base plate, for use with our Adesso Split transfer table (#10048) or a transfer table with Splitex¹ profile

10250 Adesso Multisplit transfer stand

Figure shows: *Adesso Multisplit transfer stand* #10250 with transfer table #10048 and mounted joint support



Figure shows the transfer jig #10047 in the Arto articulator with 126 mm construction height. The incisal table support is swivelled under the incisal plate

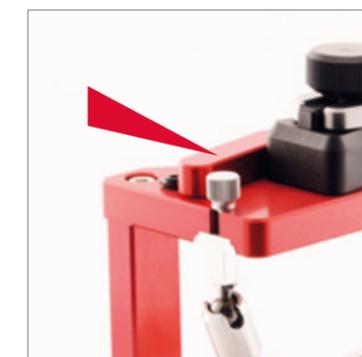


Figure shows the transfer jig #10047 in the Arto articulator with 116 mm construction height. The incisal table support is swivelled to the side

Transfer carriage for Artex¹ & Arto face bow registrations

The design enables the registration to be held directly in the articulator, without having to use the transfer stand. The rotating incisal table support is variable, making it suitable for Arto articulators 116 & 126 mm, as well as Artex¹ and Artex¹ Carbon articulators.

10047 Transfer carriages



#10049

#10050

#10250

Transfer accessories

Various Adesso Split adapters for transfer stands, positive plates, transfer tables and other accessories can be found here.

10018	Adesso Split adapter for Reference transfer stand to facilitate model transfer between the Reference and Artex ¹ systems. Adjustments are made in your Reference articulator, and the adjusted adapter is then inserted into the Reference AB transfer stand.
10043	Adesso Split positive plate suitable for subsequent conversion of a transfer table
10046	Adesso Split adapter for transfer stand SAM ² with Artex ¹ / Arto face bow for use e.g. with our universal transfer stand (#10049) to transfer the Artex ¹ face bow registration in the correct position in a SAM ² articulator.
10048	Adesso Split transfer table for transfer stand to fix a face bow registration in the articulator or transfer stand. Suitable for all articulators, transfer stands and transfer carriages fitted with Adesso Split or Adesso Multisplit.
10051	Bite fork support infinitely height-adjustable to support a face bow registration in the articulator or plaster articulator. The magnetic lower part means the bite fork support fixes itself onto all common magnetic plate systems.



Workshop service & special conversions

Our workshop service works in a *fast and affordable way*. No matter whether the work concerns coordinating or overhauling articulators, or special conversions for special applications, we will be happy to provide fast, professional assistance.

Converting scan fixators/articulators to all possible magnetic plate systems is also part of our daily work. You will get your upgraded or converted equipment back within a few days.

This is what we offer you in detail:

- Calibration of used and new articulators
- Repair service for all common articulator types
- Integration of articulator extensions
- Conversion of scan fixators/articulators
- Custom-made products and special conversions
- And much more...

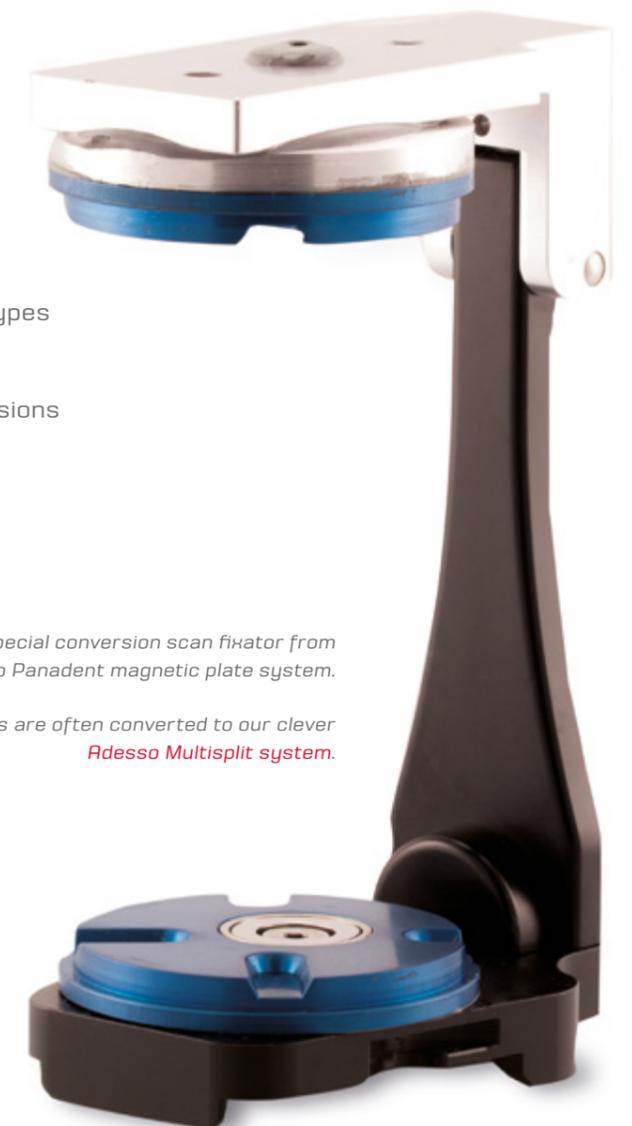


Figure shows: Special conversion scan fixator from AmannGirrbach to Panadent magnetic plate system.

These fixators are often converted to our clever *Adesso Multisplit system*.

*...do you still have questions?
We look forward to your call...
Fon +49 7236-9 33 69 0*

Sales & Service

We have limited the product information in our catalogue to a few fundamental properties and application tips. Detailed information can be obtained in the individual brochures/instruction manuals and on the Internet at: www.baumann-dental.de

Do you need a personal on-site consultation?

Give us a call: +49 7236-9 33 69 0

Pricing

Our prices can be found in our separate price list. All prices are net prices, stated in euros (€). The respective valid statutory value-added tax is stated separately. We generally deliver ex works, i.e. plus postage and packaging.

- Quantity discounts/graduated prices available on request
- We provide a 3% cash discount for prompt payment by direct debit
- All prices are subject to change. Errors and omissions excepted.

Order data

The item number always refers to 1 piece. Exceptions from this are packaging units, which are indicated accordingly.

Legal information

1 = *Artex, Splitex and Giroform* are registered trademarks of *Amann GIRRbach GmbH, 75177 Pforzheim, DE*

2 = *SAM* is a registered trademark of *SAM Präzisionstechnik GmbH, 82131 Gauting, DE*

3 = *KaVo and Protar* are registered trademarks of *Kaltenbach & Voigt GmbH, 88400 Biberach, DE*

4 = *Zeiser* is a registered trademark of *Zeiser Dentalgeräte GmbH, 71282 Hemmingen*

5 = *Quicksplit* is a registered trademark of *Hans Roßner & Sohn GmbH, 87700 Memmingen, DE*

6 = *Hipmount and Headlines* are registered trademarks of *R. Schöttli, 91052 Erlangen, DE*

7 = *Dentatus* is a registered trademark of *Dentatus International AB, Hägersten, SE*

8 = *3Shape* is a registered trademark of *3Shape A/S, DK-1060 Copenhagen K, DK*

Adesso Split, Adesso Multisplit, Model System 2000, Quadra Base,

Arto, Combiflex and Connect are registered trademarks of *K. Baumann, 75210 Keltern, DE*

General Terms and Conditions of Baumann-Dental GmbH

§ 1 Offers and placing of orders

1.1 A contract is exclusively concluded upon our written confirmation or delivery of the goods that determine the content of the contract. Insofar as the contract is explicitly confirmed by Baumann-Dental, the shipping of the goods applies as acceptance of the order.

1.2 Subsidiary agreements, amendments and additions to the contract must be confirmed in writing by us in order to be valid. The same applies to property guarantees and the waiving of the requirement for the written form.

§ 2 Terms of payment, prices

2.1 The prices are understood to be in euros ex works excluding all incidental expenses and the statutory VAT, and are stated in euros. Goods directories and prices are always subject to change.

2.2 If the payment due date is exceeded, the customer is obliged to pay interest on the purchase price at the standard bank interest rate.

2.3 In the event of worsening of the financial status and/or the liquidity of the customer following conclusion of the contract, we are entitled to exercise our rights as per section 321 of the German Civil Code (BGB) or the right to withdraw from the contract, in whole or in part, as we choose.

§ 3 Right of retention, assignment

3.1 The customer cannot refuse or withhold his services due to any counterclaims, or offset these with counterclaims, unless these counterclaims are uncontested by us in terms of their reason and extent, or are legally established.

3.2 The customer's claims against us can only be assigned with our consent.

§ 4 Delivery time, exceeding of deadline

4.1 The delivery period begins when the customer receives our order confirmation.

4.2 The delivery time shall be considered observed if, by its expiration, the deliverable has left our company or the customer has been informed of its readiness for despatch.

4.3 The delivery period is extended as appropriate in the event of measures involving industrial action, particularly strikes and lock-outs, as well as upon the occurrence of unforeseen obstacles that lie beyond our will, insofar as such obstacles verifiably have an inherent influence on the completion or delivery of the deliverable. This also applies when these circumstances occur with suppliers. We shall also not be held responsible for the aforementioned circumstances if they occur during a delay that already exists. Any contractual penalty agreed does not apply as forfeited under these circumstances.

4.4 The observance of the delivery period presupposes the fulfilment of the customer's contractual obligations.

§ 5 Passing of risk and acceptance, reporting and examination obligations, inspection

5.1 The risk passes at the latest with the shipping of the delivery parts to the customer, even if partial deliveries are made or we have assumed other services. Returns are made at the customer's risk.

5.2 Partial deliveries are permitted. § 8.1 takes priority.

5.3 The customer is to inspect the deliverable extensively immediately and report any defects, incorrect, excessive or insufficient deliveries, with proof, in writing at the latest within 10 days. The deliverable is considered approved once the period elapses. The burden of proof for hidden defects is borne by the customer.

5.4 Acceptance can only be refused due to considerable defects.

5.5 If the customer is in default with accepting the goods, we can, after setting a suitable period of notice of 14 days, withdraw from the contract or store the contractual goods at the risk and expense of the customer. Furthermore, we are entitled to auction the contractual goods following prior warning or sell the goods at the expense and for the account of the customer. The risk of accidental destruction or accidental deterioration of the contractual goods transfers to the customer at the time in which he becomes in default of acceptance.

§ 6 Retention of title

6.1 We retain possession of the deliverable until repayment of all liabilities from the business relationship (also the future business relationship) (including all outstanding balances from the current account).

6.2 The customer may undertake normal business operations and transfers to us, upon combination or mixing with other transactions, the proportionate joint ownership of the new item. Upon sale, the customer assigns to us his purchase price claim in accordance with our joint ownership share in advance as security, amounting to our invoice sum. The assignment shall be accepted by us.

6.3 In the event of conduct on the part of the customer that is in breach of the contract, in particular upon default in payment, we are entitled to withdraw the product following warning and the customer is obliged to surrender it. The exercising of retention of title and the seizure of the deliverable by us do not apply as a withdrawal from the contract, insofar as the instalment payment law does not apply.

§ 7 Liability for delivery deficiencies

For delivery deficiencies, which also include the lack of explicitly ensured properties, we are liable, to the exclusion of further claims without prejudice to section 8 as follows if the deliverable is not approved as per section 5.3.

7.1 All such parts are to be repaired or redelivered as we see fit free of charge if they verifiably display a defect that is our responsibility within 6 months from delivery and this defect has been reported to us immediately in writing. Replaced parts become our property. For significant third-party products, our liability is restricted to the assignment of the liability claims to which we are entitled against our suppliers.

7.2 No guarantee is granted for damages that are the result of the following: natural wear and tear, improper use, defective mounting or commissioning or handling, insofar as they are not attributable to our fault.

7.3 To undertake all of the repairs and replacement deliveries we consider appropriate as per our obligation, the customer is to provide us with the necessary time and opportunity with our agreement, otherwise we are relieved of liability for defects.

7.4 The warranty for the replacement part and repair amounts to 3 months, although it shall extend at least until expiration of the original warranty for the delivered item.

7.5 In the event that the customer or third party improperly undertakes amendments or repair work without our prior consent, no liability is accepted for the consequences of this.

7.6 Additional claims of the customer, particularly a claim to compensation for damages that do not occur to the deliverable itself, are excluded.

§ 8 Right of the customer to withdrawal and other liability

8.1 The customer can withdraw from the contract if it is ultimately impossible for us to render the entire performance before the passing of risk. The same applies in the event of our inability. The customer can also withdraw from the contract if, upon ordering the same type of items, the execution of part of the delivery becomes impossible in terms of quantity or he has a legitimate interest in the rejection of a partial delivery; the customer can thus correspondingly reduce the return service following mutual agreement.

8.2 If the impossibility occurs during the default in acceptance or due to the fault of the customer, the customer remains obliged to render a return service.

8.3 The customer has a right to withdraw if we allow a suitable grace period set by us for the performance of repairs or replacement deliveries with regard to a defect for which we are responsible, within the meaning of the delivery conditions, to elapse three times through our fault. The customer's right to withdraw also exists in other cases where we fail three times to provide a repair or replacement delivery.

§ 9 Copyright

9.1 Our drafts, samples, models, tools and similar are our intellectual property and may neither be imitated nor used in any other way for replication by the customer, even if no particular protective rights exist in this regard, nor are they to be given to third parties for free or for a fee, or used in any other manner than that stated as part of the agreement.

9.2 The customer is liable in the event of ordering according to the customer's requirements (including drawing, samples). In this case, trademark rights are not infringed upon.

§ 10 Place of jurisdiction, law, scope of validity for the T&Cs

10.1 German law applies. The place of fulfilment is the headquarters of Baumann-Dental GmbH.

10.2 Any conflicting terms and conditions are invalid, even if the subject matter is not governed in these terms and conditions.

§ 11 Overseas business

The following conditions also apply to overseas business:

1. All transactions are subject to German civil and commercial law. The regulations of private international law are explicitly excluded.
2. Overseas deliveries are made in exchange for payment in advance.
3. Deliveries take place in accordance with INCO Term 2010 ex works Keltern.
4. In the event of a warranty claim, service is performed following consultation with Baumann-Dental.
5. Any customs charges and bank fees shall be at the client's expense.

§ 12 Validity of the provisions

Should individual provisions be or become invalid or contain a loophole, the legal validity of the remaining provisions shall remain unaffected by this. In place of the invalid provisions, a valid provision shall be agreed that comes as close as possible to the economic purpose of the original provision; the same shall apply in the event of a loophole.



Baumann-Dental GmbH
Frankenstraße 25
75210 Kelttern / Germany

info@baumann-dental.de
www.bbaumann-dental.de

Tel +49 72 36-933 69 0
Fax +49 7236-933 69 99

Managing Directors
Karsten Baumann
Gerhard Buchenroth
Wolfgang Dawedeit

Bank details
BW-BANK
BIC: SOLADEST600
IBAN: DE35600501010004971355

Catalogue design
Alexander Kassel
(Creative Director / Baumann-Dental)